

STATMATIC

Automatic handpiece maintenance system

Worry Free™ Technology for your handpieces.



SciCan

A HIGHER STANDARD

THE ULTIMATE CARE SYSTEM FOR 3 HANDPIECES AND THEIR CHUCK MECHANISMS



Prod. no. S30000-NA

QUICK AND THOROUGH:

- > Automatically cleans and maintains up to three handpieces. Cycle time is only 15 seconds per handpiece.
- > Unique chuck care system and worry free™ leak stop feature.
- > Chemical and physical cleaning is accomplished with foaming care agent.
- > Optimal lubricant coverage throughout the handpiece.
- > Compressed-air connection allows excess lubricant to be purged.
- > Automated reproducible quality of care.
- > Easy spray can insertion.
- > Can coupling with unique safety valve.
- > Absorption pad for excess lubricant.
- > 10 years guarantee – only from SciCan.

CONCLUSION:

Improved working life of your handpieces through efficient cleaning and lubrication – saving you time and money.

TECHNICAL DATA OF STATMATIC



Dimensions (W x D x H):	Weight:	Electrical ratings:	Protection class:	Ambient temperature:	Sound level:	Relative humidity:	Compressed air:	Air consumption:	Compressed-air quality:
285 x 190 x 400 mm (11.4 x 7.6 x 16 inches)	8.6 kg 18.9 lbs	Rated voltage: 100-230 V Rated frequency: 50/60 Hz Rated power: max. 30 VA	1	5-40° C (41-104 °F)	60 (65) dB (A)	max. 85 %	4-6 bar (58-87 psi)	40 NL/min	Residual oil content: max. 0,1 mg/m ³ Residual dust: max 1µ; max. 1 mg/m ³ Residual water: max. 0.1 g/m ³ at -40 °C

ONE, TWO OR THREE. STATMATIC WORKS IT OUT.



STATMATIC is extremely straight forward: just flip down the door, snap one, two or three handpieces onto the adapters, close the door and press start. It's that simple.

STATMATIC will clean and lubricate each handpiece, one after the other, every 15 seconds. The system ensures that each instrument receives exactly the right quantity of STATMATIC lubricant. The precise dose of surfactant, oil and propellant mixture dynamically foams to actively clean all internal components of the handpiece. STATMATIC spray is expelled through a constant stream of air so that excess lubricant is removed. This also allows the water and air lines of the handpiece to be purged. As a result, your instruments are ideally prepared for sterilization. No more damaging debris or excess build up!

STATMATIC offers great flexibility. It automatically registers the number of connected handpieces, whether it's one, two or three. The STATMATIC comes with 2 E-type and 1 Multiflex* coupling. Contact your dental dealer to enquire about other coupling options.




CLEAN CHUCKS WITH CM CARE.

Secure the grip of dental burs with CM Care. The chuck maintenance tip ensures debris free chucks for safe and easy replacement of dental burs. CM-Care. It's simple. Just place the instrument head against the tip, press the button and seconds later the chuck is perfectly clean.

CM Care prolongs the useful life of your chuck, increases its functional capability and therefore, optimizes the safety of the instrument in use, for every treatment.

Prod. no. S30000

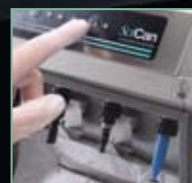

CHUCK CARE: CM CARE
Put the chuck in position,
press the button and . . . Done!

STATMATIC™ – CLEANS AND LUBRICATES STRAIGHT TO THE BOTTOM LINE.

STATMATIC is the fast and easy way to maintain your handpiece investment and extend their working life. STATMATIC makes handpiece care simple, efficient and consistent with its unique worry free™ leak stop technology.

STATMATIC is the automatic care and maintenance system that accommodates up to three handpieces and their chuck mechanisms. STATMATIC optimizes working procedures in your practice and plays an important part in cleaning, protecting and maintaining your valuable instruments. In a matter of seconds, handpieces are properly cleaned and lubricated prior to sterilization.



THE UNIQUE 10 YEAR GUARANTEE:
Highest quality stainless steel design and
engineering – you can rely on STATMATIC.
Premium Quality – Made in Germany.

You and your staff will profit from greater efficiency, fewer repairs and extended durability of your handpieces. STATMATIC pays for itself – faster than you might expect.

SIMPLE, CONVENIENT SPRAY CAN LOADING.



Prod. no. STM500-6NA

Replacing the spray can containing the maintenance lubricant is as easy as operating the STATMATIC itself. Instead of complex and cumbersome locks that need to be rotated and pushed into place, STATMATIC provides its users with simplicity and safety: just push up the retaining ring, snap the spray can into place and let the retaining ring back down again. That's it.

Thanks to its compact dimensions, STATMATIC can be located anywhere. It has its own compressed-air connection with a built-in regulator. In other words: you can connect STATMATIC directly to your compressor.

TOGETHER EVEN FASTER. STATMATIC AND STATIM.

By combining STATMATIC and STATIM™, the world's fastest autoclave, SciCan offers you the perfect solution for optimum workflow management. It's the easiest and most gentle way to clean, lubricate and sterilize your instruments.

IT COULDN'T BE FASTER:

Three instruments cleaned in 45 seconds and sterilized in 6 minutes.



CHANGING THE CAN

It's really as easy as can be: push up the retaining ring, snap the spray can into place and slide the retaining ring back down again.



CHANGING THE ABSORPTION PAD

Clean as a whistle: STATMATIC is fitted with an absorption pad to retain any excess lubricant.

Prod. no. S30102



SAFETY VALVE: WORRY FREE TECHNOLOGY

Absolutely safe: the automatic safety valve will only open when the unit has been activated.

FAX RESPONSE

If you would like more information or an appointment with a SciCan representative, simply fill out the form and fax it to SciCan. We're looking forward to hearing from you.

Yes, please send me specific information about the following SciCan products:

- | | |
|---|--|
| <input type="checkbox"/> STATIS™ Turbines | <input type="checkbox"/> STATMATIC™ Automatic Handpiece Maintenance System |
| <input type="checkbox"/> STATIS™ Handpieces and contra angles | <input type="checkbox"/> STATCARE™ Handpiece Maintenance Spray |
| <input type="checkbox"/> HYDRIM® L110w Instrument Washer | <input type="checkbox"/> STATIM® Autoclaves |

Please contact me to schedule and appointment with the regional SciCan representative

Surname, first Name

Fax

Street

E-Mail

City, ZIP

Date

Phone

Dealer Stamp

For U.S. Customers:
SciCan Inc.
701 Technology Drive
Canonsburg, PA 15317
USA
Phone: (800) 572-1211
Fax: (724) 820-1479

For Canadian Customers:
SciCan Ltd
1440 Don Mills Road,
Toronto ON M3B 3P9
CANADA
Tel.: (416) 445-1600
Fax: (416) 445-2727

